




Statement Instructions for Use for RECIPROC® Paper Points

We, VDW GmbH, located at Bayerwaldstr. 15, 81737 Munich, Germany, herewith declare the following:

For RECIPROC® Paper Points no IFU is available. Handling instructions are given on the back side of the outer packaging (example for the RECIPROC® Paper Points Assortment):

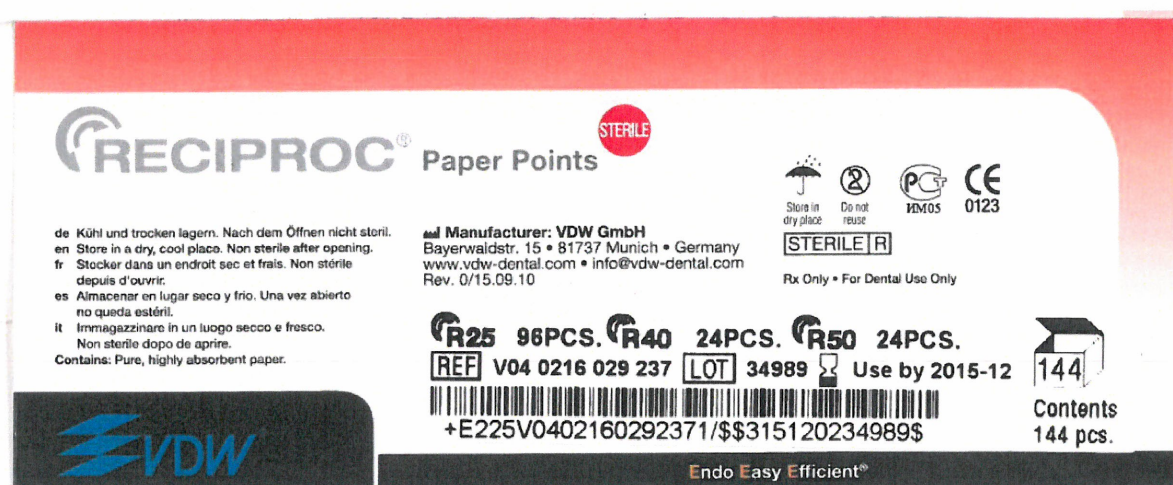

Product description

RECIPROC® Paper Points are made of absorbent paper rolled into long, narrow points with 29 mm in length. They are suitable for the RECIPROC® system. They are sterile and are available in sizes R25, R40 and R50 matching the respective RECIPROC® NiTi instrument. Length markings at 18, 20 and 22 mm help to control the working length. To support infection control the RECIPROC® Paper points are packed in blister cells of 4 paper points each that have been industrially sterilized.


Indications for Use

Paper points are used during endodontic procedures for the purpose of absorbing blood, irrigating solutions to dry root canals. They are intended only for dental, professional and single use.

Page 1 of 2

Schematic representation of the use of the RECIPROC® Paper Point:

Paper Point

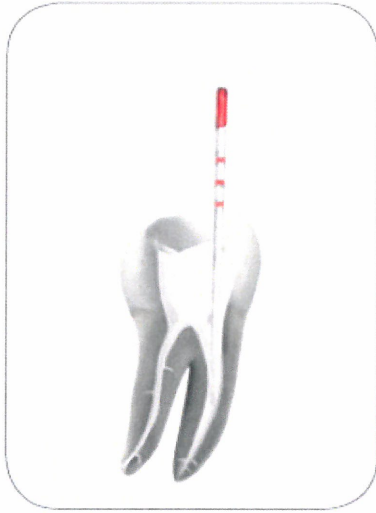


Dry the canal with the appropriately sized RECIPROC® Paper Points

Contraindications: None.

Warnings: None.

VDW GmbH
Munich, November 28, 2016

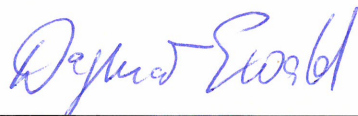

i.A. Dr. Dagmar Ewald
RA Specialist

VDW GmbH
81737 MÜNCHEN, BAYERWALDSTR. 15
TEL. 089 / 6 27 34-0, FAX 089 / 6 27 34 - 190